




Size	Waist Size
1	24" – 34" (61.0–86.4cm)
2	29" – 35" (73.7–88.9cm)
3	34" – 41" (86.4–104.1cm)
4	40" – 47" (101.6–119.4cm)
5	46" – 53" (116.8–134.6cm)
6	52" – 58" (132.0–147.3cm)
EXT.	58" + (147.32cm+)

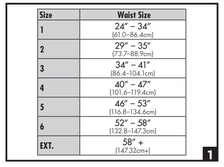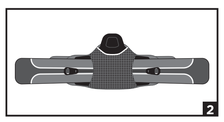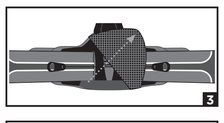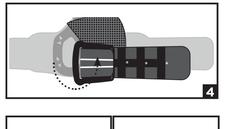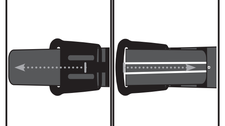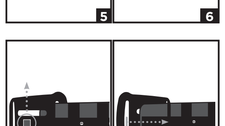
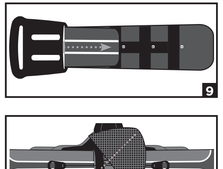
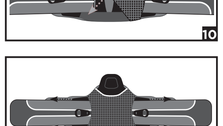
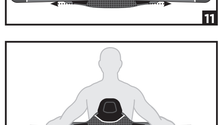
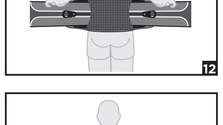
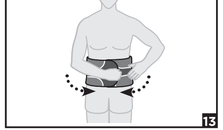
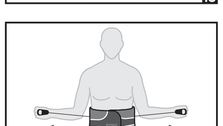
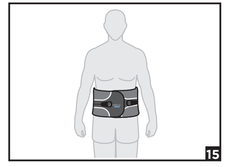
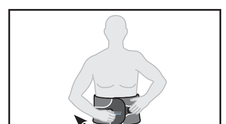
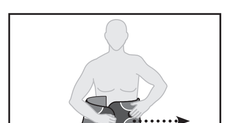
ENGLISH

INTENDED USE/INDICATIONS:

The intended use of these Lumbar products is to provide stabilization and motion restriction for the treatment of trauma related injury, acute post-surgical recovery, and the management of chronic conditions of the spine. Indicated for chronic low back pain, degenerative disc disease, severe DJD, radiculopathy, spondylosis, spinal stenosis, post-operative support, fracture management, post laminectomy syndrome, spondylolisthesis, decompression procedures, deconditioned trunk musculature, sprains/strains, bulging/herniated disc, facet syndrome, IDET procedures, SI dysfunction, sacroiliitis.

INTENDED USERS:

The intended user should be a licensed professional, the patient, or patient's caregiver. The user should be able to read, understand and be physically capable to perform all the directions, warnings and cautions provided in the information for use.

CONTRAINDICATIONS:

Not Applicable.

⚠ WARNINGS AND PRECAUTIONS:

- Carefully read fitting instructions and warnings prior to use. To ensure proper performance of the brace, follow all instructions.
- If you experience, or have an increase in pain, swelling, skin irritation, or any adverse reactions while using this product, immediately consult your licensed health care professional.
- This device will not prevent injury all injuries; these rehabilitation and activity modification are also an essential part of a safe treatment program. Consult your licensed health care professional regarding safe and appropriate activity level while wearing this device.
- Federal law restricts this device to sale by or on the order of a licensed health care practitioner.
- The frequency and duration of use shall be prescribed by a healthcare professional.
- For single patient use only.

INSTRUCTIONS FOR USE:

WAIST MEASUREMENT & SIZING

1. Measure patient's waist circumference. Use sizing chart to determine appropriate size.

WAIST BELT ADJUSTMENT

Note: Brace slips preset to Size 3 with belts in the inner slot.

1. Open brace and place it face down.
2. Lift a corner of the mesh cover.
3. Flip pulley plate over the mesh, ensuring the rigid lateral support pads remain in place.
4. Disengage belt end from pulley plate.

For Sizes 1–3:

5. Slide belt through the slot. Fold belt so its end aligns with the desired size number. Match belt to pulley plate contour, then secure to the loop.

For Sizes 4–6:

7. Remove hook piece from pulley plate.
8. Move belt piece from inner to outer slot.
9. Slide end fold the belt to align with correct size, then secure to loop.
10. Flip pulley plate back into place under the mesh flap. Tighten mesh and conceal inner slot. *Repeat steps 3 through 10 on opposite side.*
11. Position pull handles near mesh cover. Pull both pulley plates outward to center system on the back pad.

APPLYING THE BRACE

12. Align center of back pad with patient's spine and bottom edge over the lumbosacral joint.
13. Wrap right belt over left, adjusting angle for patient's body shape.
14. Tighten brace by pulling handles.
15. Attach handles to belt.

REMOVING THE BRACE:

16. Unhook pull handles and place near the back.
17. Detach front hook and loop closure and remove brace.

CUSTOMIZATION

Brace components may be trimmed, bent, or molded to fit. Refer to the product website for detailed instructions.

USE AND CARE:

Hand wash in 30°C (86°F) water with mild soap, rinse thoroughly. Rinse and lay flat to air dry. Do not machine wash or tumble dry.

MATERIAL CONTENT:

Airmesh®, Cotton, High Density Polyethylene (HDPE), Nylon, Polyester

WARRANTY:

Breg products come with a limited warranty against manufacturing and material defects. No other warranty, expressed or implied, is given and no affirmation of or by seller, by words or action, will constitute a warranty. Please see www.breg.com

DISPOSAL:

Dispose of product according to local regulations.

NOTICE TO USERS: If any serious incident has occurred in relation to the device, report to Breg and the competent authority of the Member State in which the user and/or patient is established.

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