


Size	Waist Size
1	26" – 32" (66–81.3cm)
2	31" – 37" (78.7–94cm)
3	36" – 43" (91.4–109.2cm)
4	42" – 49" (106.7–124.5cm)
5	48" – 55" (121.9–139.7cm)
6	54" – 60" (137.2–152.4cm)
EXT.	60" + (152.4cm+)

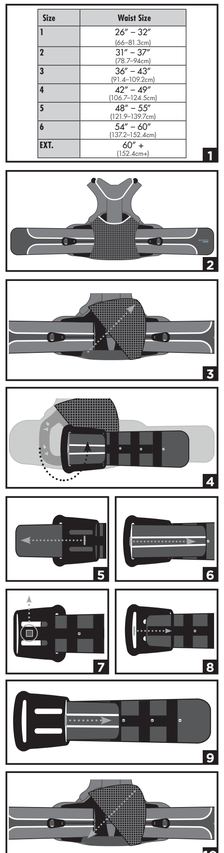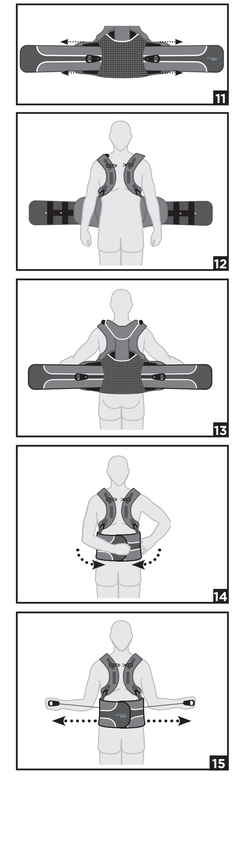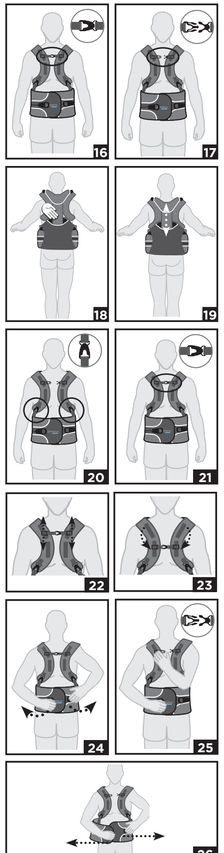ENGLISH

INTENDED USE/INDICATIONS:

The intended use of these Lumbar products is to provide stabilization and motion restriction for the treatment of severe pain related injury, acute post-surgical recovery, and the management of chronic conditions of the spine. Indicated for chronic low back pain, degenerative disc disease, severe DJD, radiculopathy, spondylosis, spinal stenosis, post-operative support, fracture management, post laminectomy syndrome, spondylolisthesis, decompressive procedures, decondtional trunk musculature, sprains/strains, bulging/herniated disc, facet syndrome, IDET procedures, kyphosis secondary to osteoporosis, scoliosis with musculoskeletal pain.

INTENDED USERS:

The intended user should be a licensed professional, the patient, or patient's caregiver. The user should be able to read, understand and be physically capable to perform all the directions, warnings and cautions provided in the information for use.

CONTRAINDICATIONS:

Not Applicable.

⚠ WARNINGS AND PRECAUTIONS:

- Carefully read fitting instructions and warnings prior to use. To ensure proper performance of the brace, follow all instructions.
- If you experience, a new or an increase in pain, swelling, skin irritation, or any adverse reactions while using this product, immediately consult your licensed health care professional.
- This device will not prevent or reduce all injuries. Proper rehabilitation and activity modification are also an essential part of a safe treatment program. Consult your licensed health care professional regarding safe and appropriate activity level while wearing this device.
- Federal law restricts this device to sale by or on the order of a licensed health care practitioner.
- The frequency and duration of use shall be prescribed by a healthcare professional.
- For single patient use only.

INSTRUCTIONS FOR USE:

MEASURING & SIZING

1. Measure the patient's waist. Use the sizing chart to choose the correct brace size.

ADJUSTING WAIST BELT SIZE

Default brace is set to size 3 with belts in the inner slot.

For sizes 1–3:

2. Open brace, lay face down.
3. Lift one corner of the mesh cover.
4. Flip the pulley plate over the mesh.
5. Detach belt end from pulley plate.
 Slide belt through slot and fold to line up with desired size number.
6. Attach belt contour to pulley plate and attach to loop.

For sizes 4–6:

7. Remove hook piece from pulley plate.
8. Remove belt panel from inner slot; insert into outer slot.
9. Slide and fold belt to align with desired size number; attach to loop.

Final Steps:

10. Flip pulley plate back into place under mesh flap.
 Reattach mesh corner, ensuring pulley plate is hidden.
 Repeat for the other side.
11. Position pull handles close to mesh cover and pull plates outward to center the pulley system.

APPLYING THE BRACE

12. Ensure pull handles are near the mesh cover and the sternum strap is unbuckled.
 Place shoulder straps on patient.
13. Align back pad with spine; position bottom edge over the lumbosacral joint.
14. Fasten waist belts (right over left), angling for body shape.
15. Pull handles away from the body to adjust tension, then secure them.
16. Buckle and tighten the sternum strap.

ADJUSTING THE POSTERIOR THORACIC EXTENSION (PTE)

17. Unbuckle sternum strap.
18. Slide hand between back pad and PTE to detach hook.
19. Adjust PTE height and reattach.
20. Ensure shoulder straps are buckled; adjust straps if needed.
21. Reattach sternum strap.
22. Slide or reposition buckles up/down on shoulder straps for height adjustment.
23. Use elastic bands to secure any extra strap length.

REMOVING THE BRACE

24. Detach the pull handles from the belt and position them close to the back of the brace.
25. Unbuckle the sternum strap.
26. Detach hook and loop at the front to remove the brace.

USE AND CARE:

Hand wash in 30°C (86°F) water with mild soap, rinse thoroughly. Rinse and lay flat to air dry. Do not machine wash or tumble dry.

MATERIAL CONTENT:

Airmesh®, Cotton, High Density Polyethylene (HDPE), Nylon, Polyester

WARRANTY:

Breg products come with a limited warranty against manufacturing and material defects. No other warranty, expressed or implied, is given and no affirmation of or by seller, by words or action, will constitute a warranty. Please see www.breg.com

DISPOSAL:

Dispose of product according to local regulations.

NOTICE TO USERS: If any serious incident has occurred in relation to the device, report to Breg and the competent authority of the Member State in which the user and/or patient is established.

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